

Dansac NovaLife FEATURE CHART

Dansac NovaLife Open
1 and 2 piece pouches with concealed closure



Dansac NovaLife Open



Dansac NovaLife 1 Open

	Mini	Mini	Midi	Midi	Maxi	Maxi
Length x width mm	220 x 132	220 x 132	239 x 144	239 x 144	259 x 158	259 x 158
Volume capacity	230 ml	230 ml	420 ml	420 ml	550 ml	550 ml
	Opaque	Clear	Opaque	Clear	Opaque	Clear
Max cut (v x h mm)	38/50 mm	38/50 mm	55/70 mm	55/70 mm	70/90 mm	70/90 mm
Starterhole	REF	REF	REF	REF	REF	REF
10 mm (cut-to-fit)	925 - 10	926 - 10	923 - 10	924 - 10	915 - 10	916 - 10
25 mm (cut-to-fit)			923 - 25	924 - 25	915 - 25	
30 mm (cut-to-fit)			923 - 30	924 - 30	915 - 30	
35 mm (cut-to-fit)			923 - 35	924 - 35	915 - 35	
40 mm (cut-to-fit)			923 - 40	924 - 40		
45 mm (cut-to-fit)			923 - 45	924 - 45		



Dansac NovaLife 2 Wafer

		36 mm	43 mm	55 mm	70 mm
Length x width mm		85 x 72	105 x 89	105 x 89	125 x 106
Ring size	Starterhole	Max Cut	REF	REF	REF
36 mm	10 mm (cut-to-fit) 20 mm (cut-to-fit) 25 mm (cut-to-fit)	28 mm	1736 - 10 1736 - 20 1736 - 25		
43 mm	10 mm (cut-to-fit) 25 mm (cut-to-fit) 30 mm (cut-to-fit)	35 mm		1743 - 10 1743 - 25 1743 - 30	
55 mm	10 mm (cut-to-fit) 35 mm (cut-to-fit) 40 mm (cut-to-fit)	47 mm		1755 - 10 1755 - 35 1755 - 40	
70 mm	10 mm (cut-to-fit) 45 mm (cut-to-fit) 50 mm (cut-to-fit)	62 mm			1770 - 10 1770 - 45 1770 - 50



Dansac NovaLife 2 Open

	Mini	Mini	Midi	Midi	Maxi	Maxi
Length x width mm	220 x 132	220 x 132	239 x 144	239 x 144	259 x 158	259 x 158
Volume capacity	280 ml	300 ml	440 ml	480 ml	650 ml	690 ml
	Opaque	Clear	Opaque	Clear	Opaque	Clear
Ring size	REF	REF	REF	REF	REF	REF
36 mm	1319 - 36	1320 - 36	1315 - 36	1316-36		
43 mm	1319 - 43	1320 - 43	1315 - 43	1316-43	1321 - 43	1322 - 43
55 mm	1319 - 55	1320 - 55	1315 - 55	1316-55	1321 - 55	1322 - 55
70 mm			1315 - 70	1316-70	1321 - 70	1322 - 70

CE
Marked products meet strict European safety and performance criteria (MDD 93/42).

Dansac NovaLife Open Features and Benefits

1 and 2 piece pouches with concealed closure

- ① **Simple and manageable opening**
The flexible closing plates enable the user to open and empty the pouch in an easy and hygienic way by applying gentle pressure on both sides.




- ② **Simple cleaning**
The staggered closing plates give easy access to the inside of the outlet for thorough cleaning if this is required.



- ③ **Simple closure**
The outlet is simply and easily folded and securely fastened by the unique hook-lock.

- ④ **Concealed closure**
The outlet is discreetly concealed into the shape of the pouch, making it more comfortable and difficult to see under clothing.

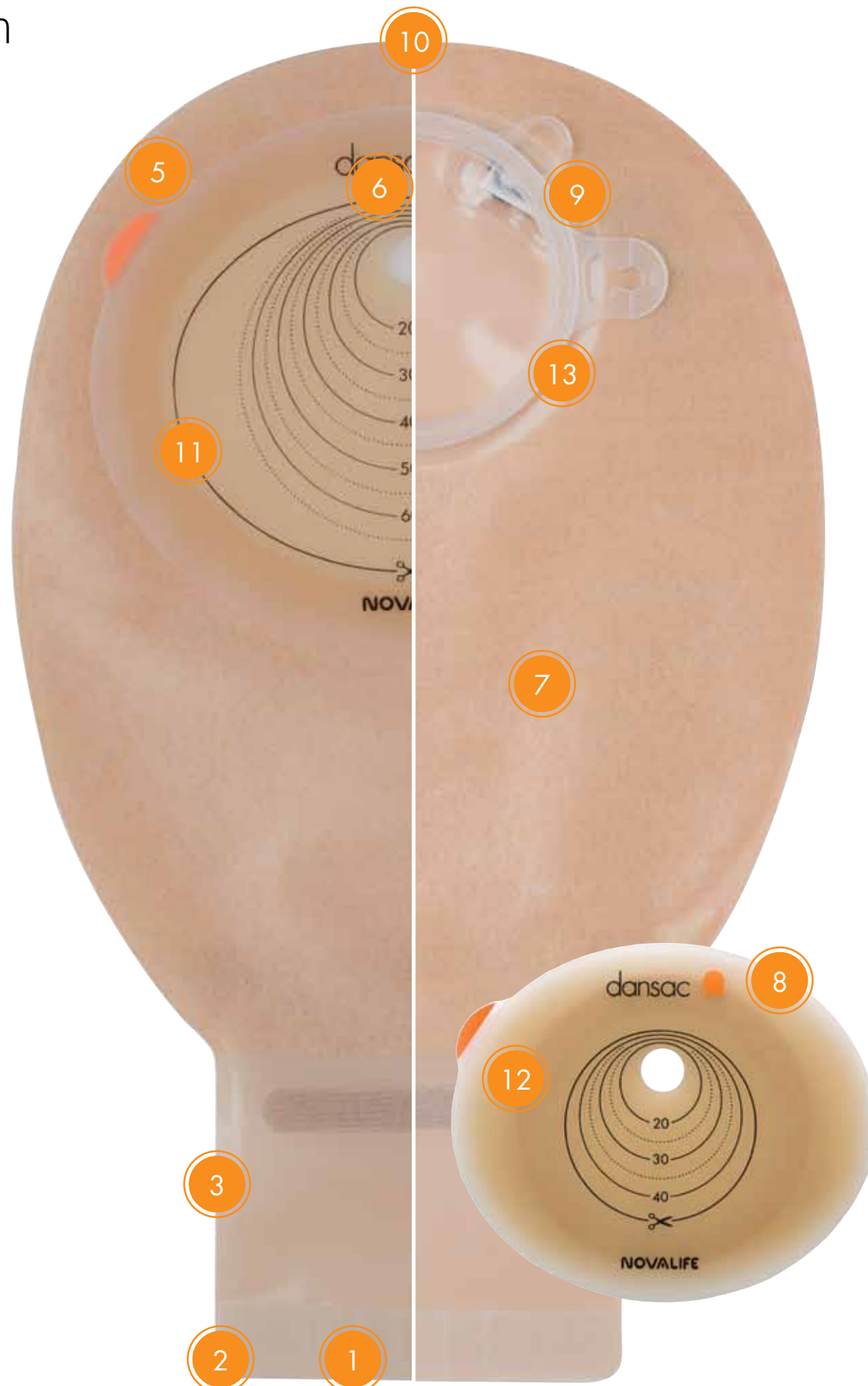

- ⑤ **Release tab**
Easy-to-see orange non-adhesive release tab.
- ⑥ **Off centre starter hole**
The pouch sits lower on the abdomen, making it easier to hide.
- ⑦ **Water repellent cover**
Stay dry even when bathing.
- ⑧ **Oval shaped barrier**
Ensures that the adhesive area is positioned where it is needed most.
- ⑨ **Highly effective Novalife filter**
The filter can expel a large amount of gas very quickly, which minimises the risk of the pouch ballooning.



- ⑩ **Low headspace**
Reduces the height of the pouch thus makes it more discreet and prevents it from slouching.
- ⑪ **Secure adhesion**
The surface of the barrier is thickest around the stoma then thins towards the edges to provide secure and gentle adhesion to the skin.



- ⑫ **Skin-friendly barrier**
GX Hydrocolloid barrier is based on natural ingredients and is proven over many years.
- ⑬ **Secure and flexible coupling**
The ring moves freely to enhance the way the wafer follows the contours of the body.